






ODS ALIGNERS

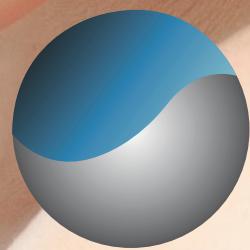

Creating Beautiful Smiles

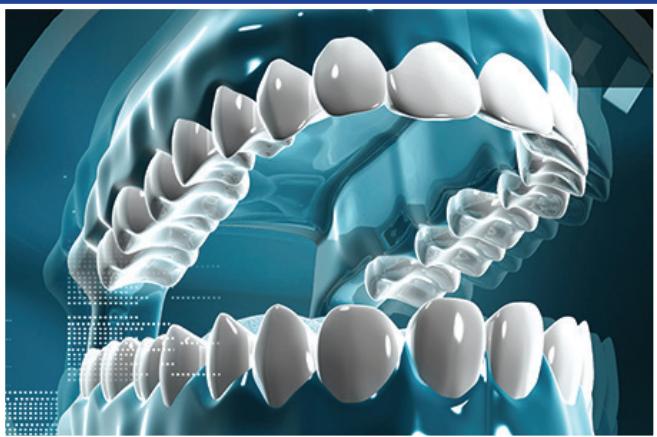

Advanced German Technology

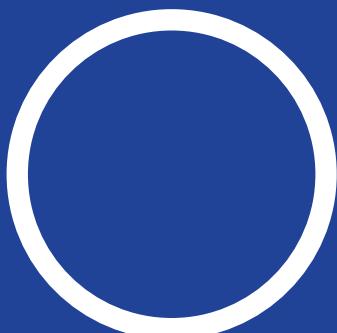

Patient-Friendly



Precise Results

# System and Advantages

**Intelligent aligner technology for visible and accurate results!**

The ODS Binary Aligners system has set new standards in orthodontic treatments. It is the fruit of years of practical experience combined with biocompatible high-performance materials. The result is an invisible, highly-effective and precise correction for almost every indication. The ODS Aligners® goes one step further with two innovations for visible and accurate results: dual-component material teamed with the Binary Aligners® technology.

In several researches and clinical trials done on various clear aligner systems, it was found that use of single material type & thickness poses a considerable limitation in orthodontic treatment done with clear aligners. One must always remember that clear aligners are, after all, orthodontic appliances and hence the principals of orthodontic treatment which applies to regular braces are applicable even to clear aligners. An important observation about orthodontic treatment is that the thickness and shape of the wire are never same. Infact it gradually increases in thickness. Similarly, the wire is never of same composition, but it changes from flexible NiTi wire to more rigid stainless Steel wire as the treatment progresses. In Clear aligner therapy, instead of metals, plastics are used and there are several types of plastics that can be used like TPU, Polycarbonate, PET, Copolyester etc as per the stage of orthodontic treatment.



# #ODSAligners

*Wide Choice Of  
Treatments*

*Accurate  
Results*


*Binary Aligners® Technology*

*Optimal Treatment  
Monitoring*

*Guarantees The  
Highest Quality  
Standards*

# What Is Our Material?

ODS Aligners use a molecularly modified copolymer plastic. Its composition transforms the inner walls of aligner into a softer and more pliable material, making it easier and more comfortable for the patient. The material is unique engineering combination of elasticity with rigidity a perfect balance.

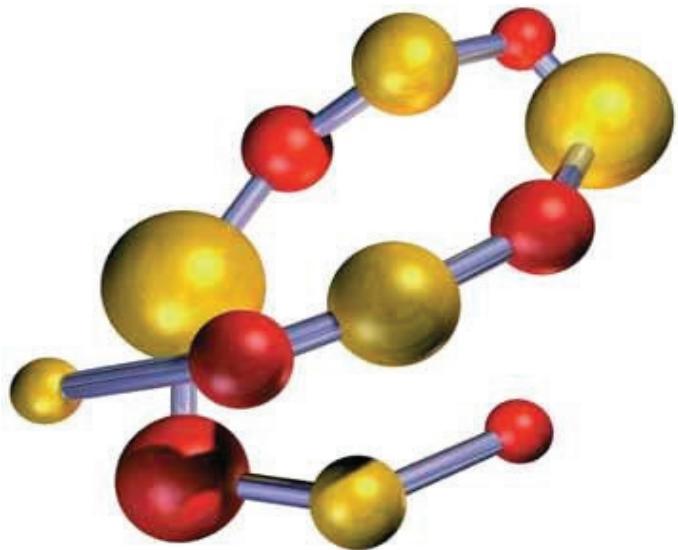

- **MAX CLARITY-** (ALL CASES) 93% light transmission, highest amongst clear aligner systems.
- **NO BREAKAGE-** as high modulus of elasticity highest in its class.
- **MORE LOAD RESISTANCE**
- **ULTRA CLEAR** final product

# *Our Offer*

**“ Win Hearts  
With Your Smile ”**



## *For Every Orthodontic Problem, ODS Has A Solution!*

From simple aesthetic correction to more complex malocclusions, ODS Aligners® have a solution for nearly 90% of cases, thanks to their system specifically tailored to suit patients' requirement. Each patient therefore benefits from an accurate and customized treatment.

Packages	Ultra	Non Complex Lite	Non Complex	Complex Lite	Complex	Complex Extended
<b>Duration (in months)</b>	Up to 3	3-6	6-10	10-15	15-24	24 & above

### **Inclusions:**

1 VDS + 1 Starter + 1 Template + All Aligners + 5 Refinement Aligners +1 Retainer

\*The 5 refinement aligners are complementary and included in the treatment duration.

# *Steps To Submit A Case*

**Step 1:** Case Submission  
(Online/Offline)



**Step 2:** Send Impressions



**Step 3:** Customized Plan



**Step 4:** Approval & Payment




**Step 5:** Aligner Production



**Step 6:** Shipping of ODS  
Aligners to the registered  
address



# *Treatment Planning*

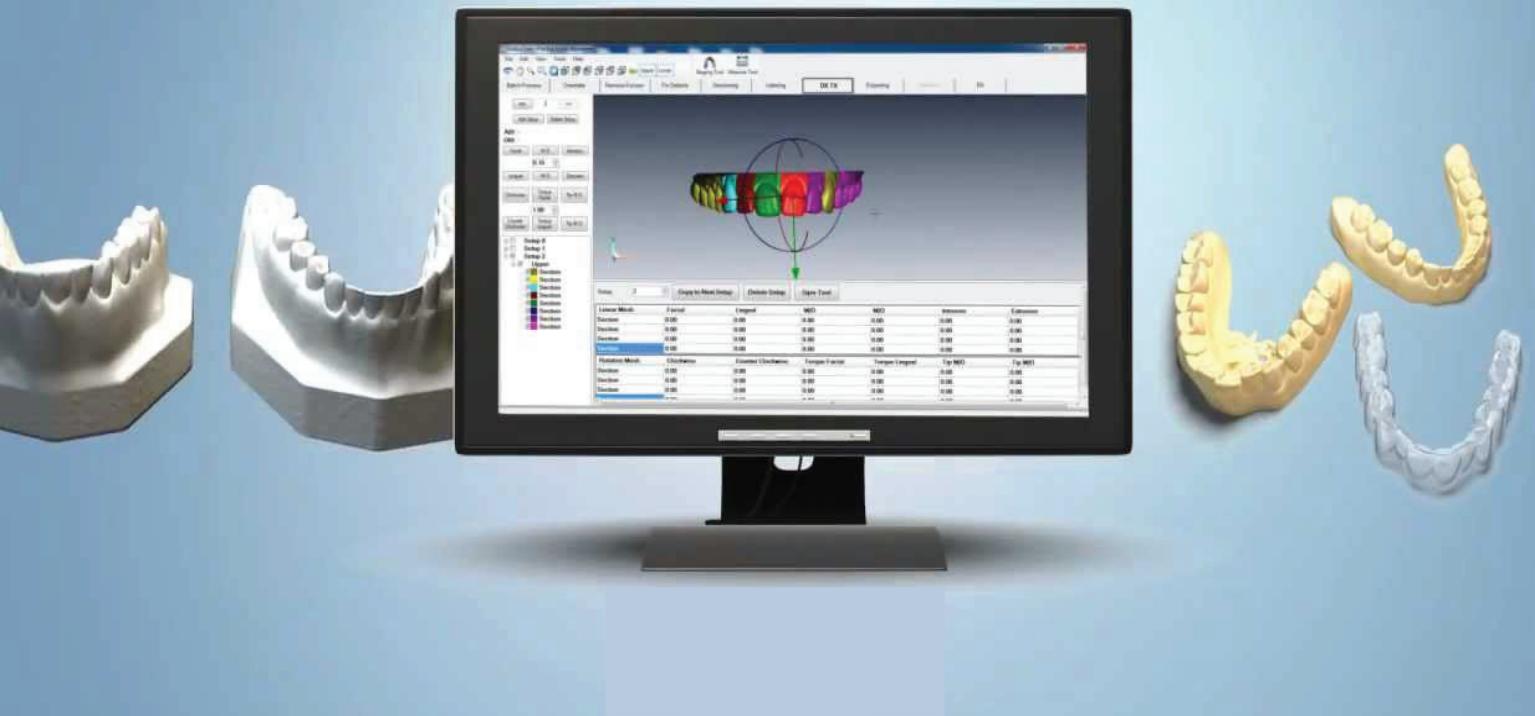

**ODS Uses Most Advanced Software For Treatment Planning**

**Simple treatment plan with practice-oriented features**

**Visualize the end result before starting the treatment**

**Explore the best possible treatment option**

**Helps to ensure every required tooth movement**

**Zero chance of any error in treatment planning**

# Patient Record

For a case, we require:

- Impression/Digital Scan
- Intra & Extra Oral Images
- X-Rays
- Dental/Medical History
- Patient Local/Travelling



## Intraoral Images

Proper Images



Improper Images

# Attachments/Engagers

## What Are Attachments/Engagers?

Aligner attachments are tooth colour small blocks used to provide anchorage, movements.



*Attachments  
(without aligner)*



*Attachments  
(with aligner)*

## Use of Attachments/Engagers

- Provide the optimal force system to achieve predictable tooth movements.
- Are designed for the width, long axis, and contour of each tooth.
- Are positioned precisely to deliver the forces while eliminating interferences.
- Crown size and tooth anatomy will determine the ideal placement.
- Are used individually or in predefined combinations.

# Attachments/Engagers

## BE CONFIDENT IN YOUR SMILE!

LET US HELP YOU ACHIEVE &  
MAINTAIN YOUR BEAUTIFUL SMILE!

### *Types of Attachments/Engagers*

Shape	Feature	Available On	Movement
	Conventional	All Teeth	Retention
	Rotational Dimples	Incisors & Canines	Couple force For Rotation
	Buccal Torque Ridge	Incisors	Lingual Root Torque
	Buccal+ Lingual Torque Ridge	Upper Incisors	Lingual Root Torque+ Retraction
	Optimized Root Up righting	Canines & Premolars	Root Up righting
	Optimized Rotation	Canines & Premolars	Rotation
	Optimized Anchorage	Maxillary Molars & Second Premolars	Anterior Retraction
	Horizontal Bevelled	Incisors	Extrusion

To know more about attachments visit: [www.odsaligners.com](http://www.odsaligners.com)

# *Successful Cases*

## *Category: Ultra*



*Before*



*After*



*Before*



*After*



*Before*



*After*

# *Successful Cases*

## *Category: Non-Complex*



*Before*



*After*



*Before*



*After*



*Before*



*After*

# *Successful Cases*

## *Category: Complex*



*Before*



*After*



*Before*



*After*



*Before*



*After*



**CONSULTING**



**SOLUTIONS**



**CONFIDENCE**



**ASSISTANCE**



**INFORMATION**



**INSTRUCTIONS**



**TIPS**



**SERVICE**

### **Advice and Support From Our Team of Experts**

At ODS Aligners®, we are committed to our client- practitioners and are devoted to providing you with solutions that really matter– every day in your practice and clinic, you can count on our support. Our network of experts (orthodontists, dental technicians and consultants) is here to support you.

### **Supporting Lines**




**Telephone**



**In Person**




**Online**

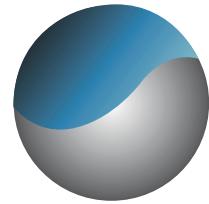

## ODS Presence Across India

- Head Office
- Providers
- Franchisees
- Upcoming Franchisees

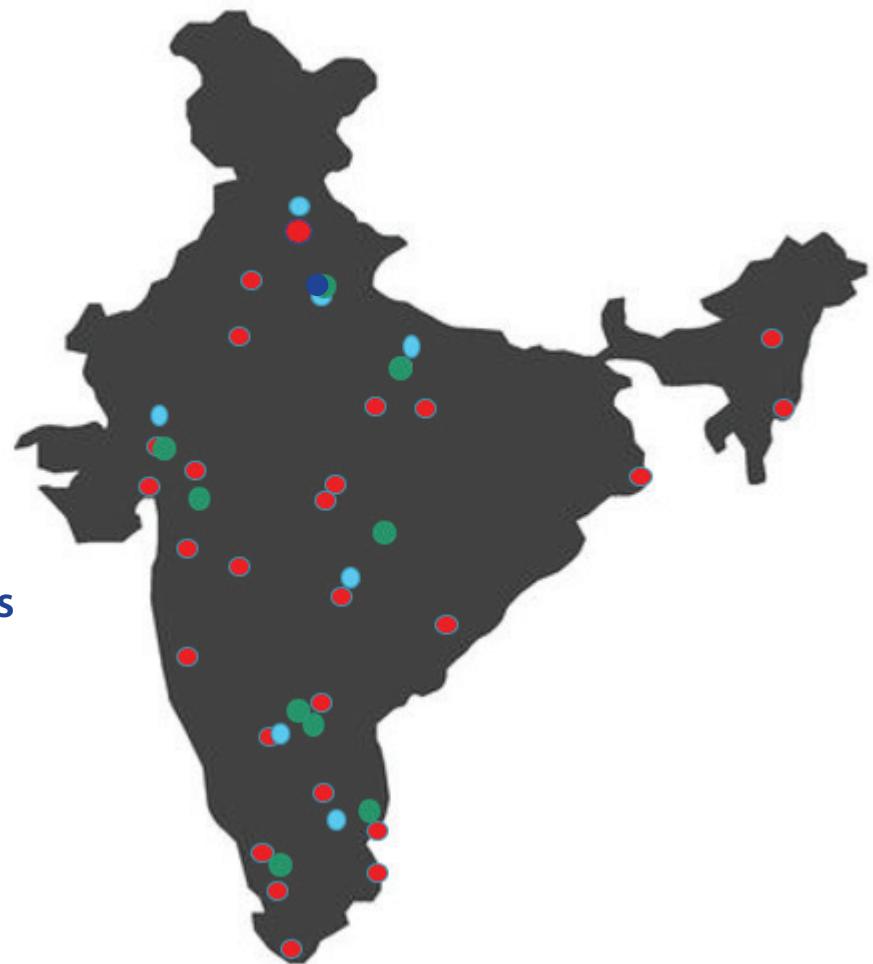

**Orthodontic Solutions Pvt Ltd**  
(Ph)+91-755 755 2921, 87440 13222/666  
Landline: 0120-4164447



/odsaligners



/odsaligners



/odsaligners



/odsaligners